

Trouble hearing?

Do you have trouble hearing in some situations?



At home. Is it hard to follow the conversation when you are chatting with a friend or relative?



At work. Do you have trouble hearing what's being said at meetings or seminars?



At school. Do you find it difficult to hear what your teacher is saying in class?



Outdoors. Is it hard to hear what your partner is saying when you are outdoors for a walk?



Watching TV. Do you miss words when watching TV at normal volume?



With hearing aid. Do you feel that your hearing device doesn't fully help in some situations?

This is the solution

Superior sound. Domino is a revolutionary personal hearing system that uses digital state-of-the-art sound processing and wireless wideband technology to bring out speech in difficult listening situations. Domino delivers crystal-clear stereo sound, free from background noise and feedback, and with no risk of eavesdropping.

Whether you only have trouble hearing from time to time or are already using a hearing device, it will make a big difference in your daily life.

Easy to use. Domino's unique design makes use very intuitive. All common features can be accessed with a single push of a button, eliminating complicated menus.

High quality. Domino is designed by Bellman & Symfon of Sweden and comes with a two-year worldwide warranty.

Domino Classic



Discover the difference.

Discover how Domino can improve your hearing in daily life.



Packed with useful features



High-Definition Audio

Domino uses all-digital technology to stream audio in full stereo. The result is crystal-clear stereo sound with outstanding clarity and definition.



Wireless Secure Encryption

Domino uses state-of-the-art digital encryption technology, resulting in interference-free speech, absolutely impossible to eavesdrop.



Easy Connect

Domino is ready for action with the touch of a button, and you can just as easily connect it to your TV, stereo system, computer or MP3 player.



Travel Smart 2.4 GHz

Domino uses the 2.4 GHz global frequency standard. This means you can use your Domino system anywhere in the world.


How Domino Classic works

The Domino Classic system consists of a transmitter and a receiver, equipped with state-of-the-art microphones.

How it works. Sound is picked up by the transmitter and receiver microphones and converted into digital data. Speech is clarified and annoying background noise is eliminated. The sound is then streamed over an encrypted wireless radio link from the transmitter to the receiver. The transmitter also incorporates a stereo input for streaming high quality audio from your TV, MP3 player or stereo system.

Maximum control. The receiver amplifies sound and sends it to your headphones, earphones or neck loop. From the receiver, you can easily control the functions of the two units.



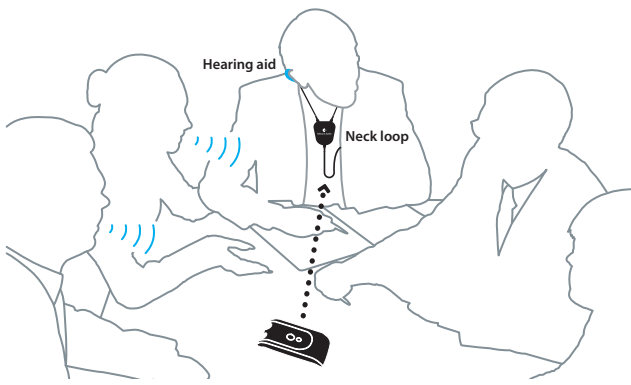

At work



Are there sometimes intense discussions around the conference table? When more than one person speaks at the same time it can be real be hard to distinguish what is being said.

Domino helps you to hear better in these situations and you can easily put the transmitter close to a distant speaker of your choice.

In school

Are you sitting in the back of the classroom? Even if your teacher is far away, Domino allows you to experience the teacher's voice as if he was sitting next to you, completely clear and noise free.

Want to talk to a schoolmate beside you? Simply push the Mic-button to hear your schoolmate's voice along with your teacher's.

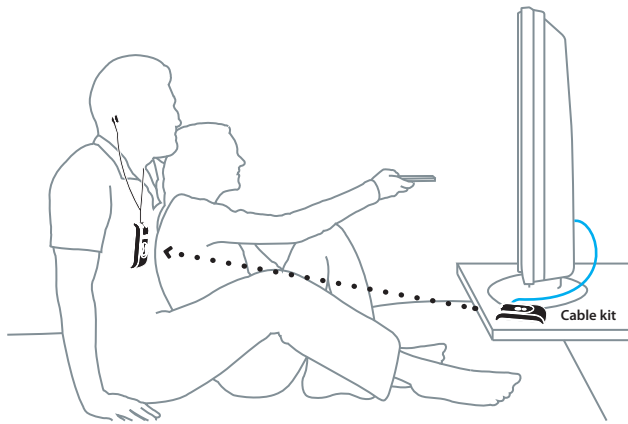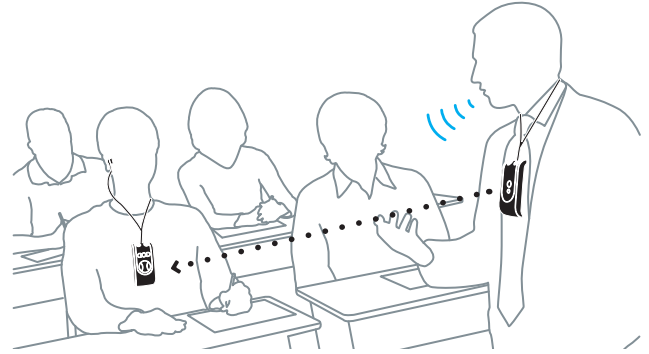

At home


Connect your Domino transmitter to the TV and sit back and enjoy clean and crisp sound in full stereo. With Domino, you don't have to crank up the volume just to follow your favorite show. Great news for your spouse and neighbors too.

When you want to discuss the film with your partner, you can easily mix the sound from the TV with the internal microphones on your receiver.

Outdoors

When you are outside taking a walk speech is often directed away from you, and there is ever so often some distance between you and your partner. This in combination with traffic noise can make it very difficult to follow a conversation.


Domino suppresses annoying background noise and overcomes the distance from the speaker very effectively.

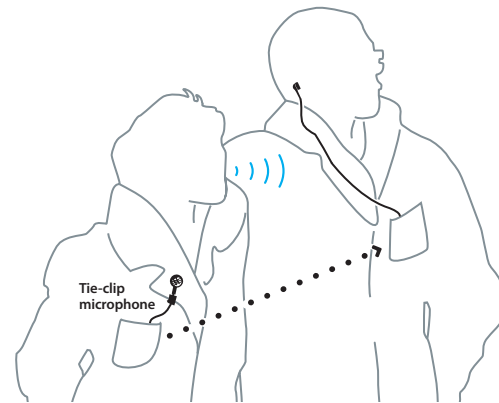

With hearing aid

Domino radically improves the performance of your hearing device, especially when you need to hear sounds at some distance. By connecting a neck loop to your Domino and switch your hearing aid to T-mode, you can take full advantage of Domino's vast advantages, for instance when attending meetings or lectures.

You can also use Domino with your hearing device when you watch TV or listen to music, very comfortable and completely wireless.

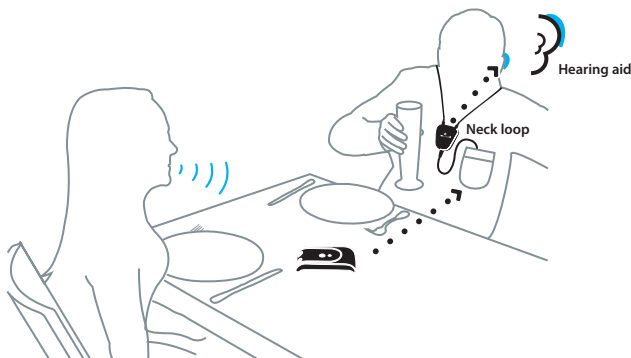

Tailor made.

Create a solution completely tailored to your needs.

Create your own solution

The Domino Classic package includes many useful accessories. If you want to create a solution completely tailored to your needs, you can also pick your choice from a wide range of original accessories. Chose Domino Classic and discover the easy way to interact with your friends, colleagues and classmates.

In the box

- Transmitter unit
- Receiver unit
- Protective travel case
- Earphones (in-ear)
- User manual
- USB universal charger with mains for EU, USA, AU & UK
- USB charger cable
- Neck strap with safety clasp (2)
- Also available with E9128 Neckloop (Art. no. BE8017)

▼ Travel case that protects the Domino units and accessories.




Original accessories



Headphones BE9122
Robust stereo headphones with great sound quality. Features semi closed drivers that reduce sound leakage. Foldable for easy travel.



Stetoclips BE9125
The Stetoclips are capable of delivering sound pressure up to 140 dB. They are comfortable to wear and very easy to put on and take off.

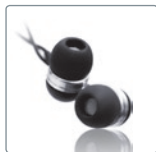

Earphones BE9124
Small and light weight stereo in-ear earphones with great sensitivity and excellent isolation from ambient noise. (Included in the package)



Neckloop BE9128
With a neck loop, Domino becomes best friend with your hearing aid. Just plug it into your Domino and switch your hearing aid(s) to T-position.



Cable kit BE9137
Use the cable kit to connect Domino to your TV or multimedia equipment and enjoy high quality stereo sound. Tele plug and RCA adapter is included.



Tie-clip microphone BE9136
Attach it close to the sound source in extremely difficult listening environments. Also a great choice for a discreet look as the transmitter can be hidden in a pocket.



USB car charger BE9143
With the car charger you can easily charge your Domino units in your car. Simply plug it in, and your Domino will be fully charged within 2.5 hours.



Extra receiver BE2250
It is easy to add one or more extra receivers to your Domino system, for example in a conference hall with multiple listeners.



Scart switch BE1360
The Scart switch automatically switch source on your various TV-equipment and transfers the signal in stereo to your Domino or hearing aid (via a neck loop).



Bellman & Symfon®

Head Office
Bellman & Symfon Europe AB
Södra Långebergsgatan 30
421 32 Västra Frölunda
Sweden
Phone +46 31 68 28 20
E-mail info@bellman.se
www.bellman.se

China Office
Bellman & Symfon Asia Ltd.
115-116, Main Tower, GD Int'l Hotel
339 Huanshi Dong Rd
Guangzhou, China, 510095
Phone +86 20 8318 2950
E-mail info.cn@bellman.com
www.bellman.cn

Contact your closest Bellman & Symfon reseller: